To Clean and Maintain the Product:

1. Cleaning

Before you clean, please make sure that the product is switched off and the power plug has been pulled out from the socket. Please use dry cloth to clean the product.

ATTENTION

Please make sure that you do not use Benzene, thinner and other aggressive detergents to clean this product. The product may fade as a result of doing so.

2. Maintenance

Make sure that the switch is off and the power plug has been pulled out of the socket. Keep the product clean and away from moisture and water.

Technology Parameter:

Name: 3 in 1 Heated Massager Rated Voltage: DC12V,2.0A

Rated Power: 24W SKU: SBA-00080 Factory Code: 828-VB



Correct Disposal of this product:

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for envirnmental safe recycling.

3 IN 1 HEATED MASSAGER USER MANUAL



This manual should be kept on hand for easy reference. It is your guide to a safe and efficient operation.

Simply Brands
www.simply-brands.com

Importer Name: Simply Brands (Asia) Ltd

Lockhart Road, Wanchai, Hong Kong

Important Safety Instructions:

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER - To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1. The product should never be left unattended when plugged in. Unplug from outlet when it is not in use and before putting on or taking off parts.
- 2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3. Close supervision is necessary when this product is used by, or near children, individuals, or disabled persons.
- 4. Use this product only for its intended use as described in this manual. Do not use attachments not recommended buy the manufacturer.
- 5. Never operate this product if it is dropped, has a damaged cord or plug, a malfunction, or it is dropped in water. Return the product to service centre for examination and repair.
- 6. Do not carry this product by its supply cord or use cord as a handle.
- 7. Keep the cord away from heated surfaces.
- 8. Never operate the product with the air openings blocked. Keep the air openings free or lint, hair and the like.
- 9. Never drop or insert any object into any openings.
- 10. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 11. To disconnect, switch off all controls and remove plug from outlet.
- 12. The product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or with the lack of experience and knowledge-unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
- 13. Children should be supervised to ensure that they do not misuse the product.
- 14. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

Save these instructions - this product is intended for household use only.

How To Operate:




- Power
- Direction
- Speed
- Heat

- Press the "Power" button, the product will start kneading and heating automatically.
- 2. Press the "Heat" button to turn on/off.
- 3. Press the "Direction" button to change the massage direction.
- 4. Press the "Speed" button to adjust the massage speed.
- 5. Press the "Power" button again to turn off the massager

Massage Positions:




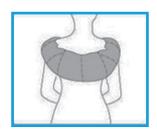


Neck

Shoulders

Upper Back

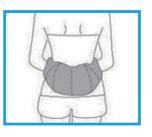

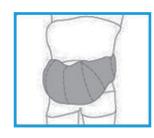



Hips

Abdomen

Thighs

Name and Function of the Machine Parts:

